

Q500[®] SERIES

Make your own adventure



Q Series
Powerchairs

Advanced Engineering. Customised for you.



Q500 R



Q500 M



Q500 F

A man with a beard, wearing a grey long-sleeved shirt and brown trousers, is seated in a black motorized wheelchair. He is holding a smartphone in his right hand and looking out over a large body of water. The sun is low on the horizon to the left, creating a warm, golden glow and reflecting on the water's surface. The background features a line of trees and hills under a clear blue sky. The wheelchair has large black wheels with orange accents and a smaller front wheel.

ADVANCED ENGINEERING. CUSTOMISED FOR YOU.

QUICKIE, JAY, WHITMYER and SWITCH-IT bring you the most comprehensive range of premium, innovative and technically advanced powerchairs, each created individually for you whatever your clinical or lifestyle need.

THE POWER OF QUICKIE®

From basic mobility to complex rehab. A wide range of bases, drive tractions, seating systems, electronics, controls and accessories, that can be combined to create the perfect powerchair for each individual need.

CLINICAL NEED



Q100

SEATING
SLING SEAT | **BASES**
RWD



Q200

SEATING
SLING SEAT | **BASES**
RWD



Q500

SEATING
SEDEO PRO | **BASES**
MWD
FWD
RWD



Q700

SEATING
SEDEO ERGO | **BASES**
MWD
FWD
RWD



Q700 UP

SEATING
SEDEO ERGO UP | **BASES**
MWD
FWD

■ RWD: Rear-Wheel Drive
 ■ MWD: Mid-Wheel Drive
 ■ FWD: Front-Wheel Drive

BASE PERFORMANCE

SERIES OVERVIEW

Q500 SERIES OVERVIEW

Resulting from years of experience of market leading drive bases, seating and advanced electronics, QUICKIE® introduces the QUICKIE Q500 Series. Explore the outdoors with the Q500 and make your own adventure!



**Q500 WITH
SEDEO PRO**
REAR-WHEEL DRIVE



**Q500 WITH
SEDEO PRO**
MID-WHEEL DRIVE



**Q500 WITH
SEDEO PRO**
FRONT-WHEEL DRIVE



ROBUST OUTDOOR POWERCHAIRS

With a choice of three drive bases (all of them with advanced suspension technologies and powerful motors), the highly customizable SEDEO[®] PRO seating (including JAY and WHITMYER) and a wide range of high end R-Net and SWITCH-IT controls, you can tailor your Q500 individually to you whatever your clinical and lifestyle needs are.



PERFORMANCE BASES

Q500 M: Mid-wheel drive

Q500 R: Rear-wheel drive

Q500 F: Front-wheel drive

SEATING SYSTEMS

SEDEO PRO

JAY & WHITMYER

CONTROLS & CONNECTIVITY

VR2 control

R-NET controls

SWITCH-IT controls

LINK-IT mounting system

DISTINGUISHED STYLE

STYLISH DESIGN and choice of colours

CHOOSE YOUR BASE

FRONT, REAR OR MID-WHEEL DRIVE

All the Q500 bases offer outstanding outdoor performance combining powerful 4 torque motors, drive wheel suspension for great comfort and traction, large wheels, speeds up to 13 kph and batteries up to 80 Ah.



The choice of a drive base affects how the chair performs in different situations and environments. QUICKIE Q500 offers the most complete range of bases in the market so each user can select the best option based on his/her personal preferences, experience, clinical, lifestyle and environmental needs.

	Q500 R	Q500 M	Q500 F
Driving type	Rear-wheel drive	Mid-wheel drive	Front-wheel drive
Turning radius	950 mm	560-630 mm	735 mm
Drive wheels	14"	13"	14"
Castor wheels	9" & 10"	7"	9" & 10"
Max. speed	6 & 8 kph (standard), 10 & 13 kph (optional)		
Battery size	60 Ah; 80 Ah optional		
Max. reversing width	1258 mm	1100 mm	1240 mm
Max. safe slope	10°	8°	8°
Max. kerb climbing	100 mm with KC 80 mm without	100 mm	70 mm
colour inlay choices	7 colours	7 colours	7 colours
Max. range (ISO 7176-4)	40 km	40 km	40 km
Lower seat to floor height	415 mm with Tilt / 425 mm with Lift & Tilt (with 60 Ah batteries)		

Q500 R
REAR-WHEEL DRIVE

The rear-wheel drive is a stable base with good kerb climbing ability.

Easy to drive, provides good control with non proportional input devices.



Q500 M
MID-WHEEL DRIVE

Mid-wheel drives offer unmatched maneuverability, with a smaller foot print than rear or front, it offers the most compact turning circle and is extremely intuitive to drive and control.



Q500 F
FRONT-WHEEL DRIVE

The front-wheel drive base offers great outdoor performance with good indoors maneuverability, and best stability on steep declines.



PERFORMANCE BASES

Q500 M

MID-WHEEL DRIVE BASE

The new QUICKIE Q500 M base balances superb indoor maneuverability with outstanding outdoor performance. Its powerful motors combined with 6 wheels suspension and QUICKIE's anti-pitch technology, will ensure you have the most comfortable and safe driving experience.

SUPERIOR INDOOR MANEUVERABILITY

Great access through narrow doorways, in lifts and around confined spaces, with the compact turning radius and base dimensions. Easily transfer into and out, effortlessly access vehicles and under tables with the lower than ever before seat to floor height from 41.5 cm with powered COG Tilt, or 42.5 cm with Lift & Tilt.



SUPERB OUTDOOR PERFORMANCE



ALL WHEELS SUSPENSION

Enjoy a smooth and stable ride, even over rough terrain, with QUICKIE's patented all-wheel suspension - six independently suspended wheels absorb all those bumps. Climb kerbs of up to 100 mm and have the best traction even when ascending gradients of up to 8 degrees.



ANTI-PITCH TECHNOLOGY

QUICKIE's anti-pitch technology eliminates pitching on slopes, and with all your weight located directly above the drive wheel to maximise traction you've never been in a safer driving position.



MAXIMUM FREEDOM

With ultra-reliable high torque 4-pole motors powering speeds of 6, 8, 10 or 13 kph and a choice of 60 or 80 Ah batteries to maximize your range, you can always travel in confidence with the Q500 M.

PERFORMANCE BASES

Q500 F

FRONT-WHEEL DRIVE BASE

Q500 F offers a lower seat to floor height of 415 mm to provide good access to wheelchair adapted vehicles, tables, and many more.

GOOD INDOOR MANEUVERABILITY

Q500 F, with in-rigged castor arms design, is responsive and offers good maneuverability for indoor use.



OUTSTANDING OUTDOOR PERFORMANCE

The Q500 F offers a robust design combined with advanced drive wheel suspension, leveller technology and powerful motors to ensure a smooth and stable ride in any terrain.



ADVANCED DRIVE WHEEL SUSPENSION

With the Full Track Suspension System you can overcome any obstacle and barely feel it. Together with the leveller technology, it will ensure you a smooth and comfortable ride in any terrain!




LEVELLER TECHNOLOGY

Maximum traction thanks to the unique leveller technology. Both drive wheels are connected by a pivot bar that absorbs the natural pull to one side when travelling across cambers and round corners at speed so both wheels keep maximum traction at any time.



MAXIMUM FREEDOM

With large wheels, powerful high torque 4-pole motors powering speeds of 6, 8, 10 or 13 kph and a choice of 60 or 80 Ah batteries to maximize your range, you can always travel in confidence with the Q500 F.

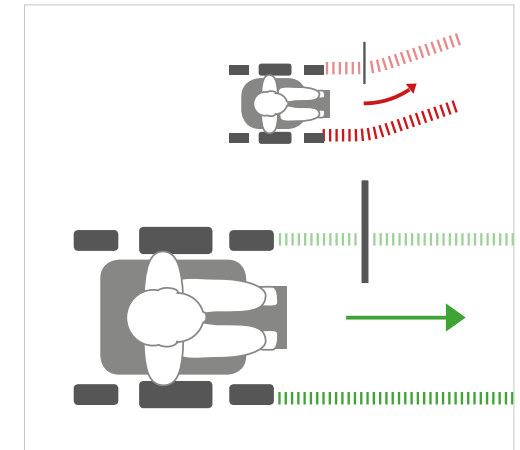

KEEP ON TRACK WITH GYRO-CONTROL

The optional Gyro-Tracking System automatically detects and corrects the smallest deviation from your direction. By providing a pin-point accurate drive, it'll reduce the number of physical or cognitive actions required to keep you on-track – making for a less-demanding drive (especially if using special controls).

PERFORMANCE BASES

Q500 R

REAR-WHEEL DRIVE BASE

Q500 R offers a lower seat to floor height of 415 mm you can access your wheelchair chair adapted vehicle too and many more

EASY AND INTUITIVE

Q500 R delivers an easy to handle and intuitive drive due to the drive wheel position.



BUILT FOR OUTDOOR PERFORMANCE



FULL TRACK SUSPENSION

With the Full Track Suspension System you can overcome any obstacle and barely feel it. Together with the leveller technology, it will ensure you a smooth and comfortable ride in any terrain!



LEVELLER TECHNOLOGY


Maximum traction thanks to the unique leveller technology. Both drive wheels are connected by a pivot bar that absorbs the natural pull to one side when travelling across cambers and round corners at speed so both wheels keep maximum traction at any time

The Q500 R Robust design combined with the unparalleled performance of Full Track Suspension System, leveller technology and built in traction control you don't need to compromise on a smooth, stable ride with full control of your chair at all times.



MAXIMUM FREEDOM

With large wheels, powerful high torque 4-pole motors powering speeds of 6, 8, 10 or 13 kph and a choice of 60 or 80 Ah batteries to maximize your range, you can always travel in confidence with the Q500 R.

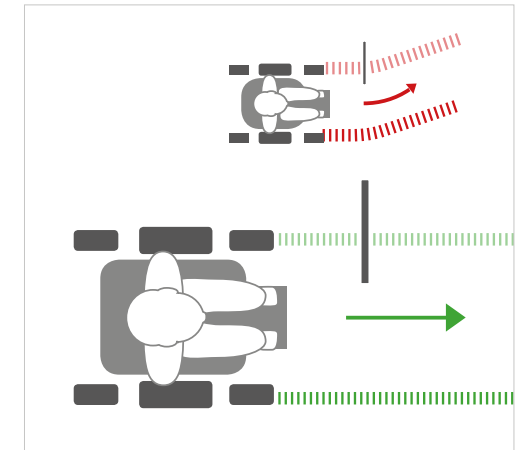

KEEP ON TRACK WITH GYRO-CONTROL

The optional Gyro-Tracking System automatically detects and corrects the smallest deviation from your direction. By providing a pin-point accurate drive, it'll reduce the number of physical or cognitive actions required to keep you on-track – making for a less-demanding drive (especially if using special controls).

SEDEO PRO

The SEDEO PRO seating system combines a highly adjustable and configurable seating frame with an excellent choice of support and positioning options, including a wide range of SEDEO PRO seat cushions and backrests for mild to medium clinical needs.

The SEDEO PRO is very comfortable to use and enhances user independence in everyday life.

QUICK & EASY TO ADJUST

The SEDEO PRO frame offers great flexibility as it is fully adjustable . All adjustments can be done in a very quick and easy way using just one tool, even with the user in the chair. It can be easily customised to suit individual requirements and quickly adapted to changing needs.



HIGHLY CUSTOMISABLE SEATING



Anti-shear backrest

Discover the new powered recline backrest with 100 mm anti-shear - makes sure that there is no sliding between you and the cushion. Additionally the SEDEO PRO is available with a fixed backrest or two other reclining backrest options: Mechanical and with actuators.

Power elevating legrest

The power elevating and length compensating legrest allows you to reach the best position. It is available as a centre mount and a swing away. Additionally you can also choose from a standard centre mount or 70° swing away legrest and a manual elevating swing away legrest.

Flexible back cushion

A wide range of wheelchair seat cushions and backrest cushions is available for the SEDEO PRO seating system, including the flexible back cushion with adjustable straps and lateral support - for unique comfort and easy positioning.

Positioning elements

Many positioning elements are available like SEDEO headrests & WHITMYER headrests, swing away lateral thoracic supports, thigh supports or knee supports. All positioning elements are easily swing away or detachable.

Flip back armrests

The easy to operate flip back armrest gives you free transfer space. You can choose from a removable or reclining version. Also available: Single post removable armrest. All compatible with 5 armpad options.

SEATING SYSTEMS



JAY[®] & WHITMYER[®]

POSITIONING OPTIONS

A wide range of WHITMYER headrests and JAY cushions are available and compatible with the SEDEO PRO frame to meet higher clinical needs and ensure we offer the best seating possible.



WHITMYER HEADRESTS

Created to address a range of user needs from basic support to multiple-point control. WHITMYER headrests come in a variety of pad shapes and sizes for superior head positioning and comfort.

JAY CUSHIONS & BACKS

Designed for pelvic stability, skin protection and comfort, JAY cushions and backs are low maintenance, highly customisable and are available in a large variety to meet your needs.



CONTROLS & CONNECTIVITY

CONTROL AT YOUR FINGERTIPS

Be the master of your environment with the Q500 Series and its advanced controller technology. Not only can you control your Q500 the way you want, but also your household appliances, phones and computers.



VR2 & R-NET

ADVANCED JOYSTICK

Every Q500 comes with a choice of the reliable VR2 control and the superior R-Net advanced joystick. A large, easy-to-read HD display and simple to control paddle switches allow for intuitive operation.

CONNECT

Seamlessly take control of your environment using the integrated Infrared functionality. Control your home entertainment system, TVs, DVDs, docking stations and other home entertainment systems simply through your joystick.

BLUETOOTH & IDRIVE

Control your phone, computer, tablet and any other Bluetooth and i-devices with this useful optional extra. Connect and control up to four devices.

QUICKIE'S PATENTED ASSIGNABLE BUTTONS

Change the function on any button on your joystick to any command you desire, thanks to QUICKIE's patented assignable buttons. You can also assign a second (long press) function to buttons, giving you quicker, simpler and smoother access to functions you commonly use without having to use your joystick.



SWITCH-IT[®]

INNOVATIVE CONTROLS

Conveniently access your 5 programmable positions with the Ctrl+5. It allows you to go straight to the function you need at the touch of a button.

ALTERNATIVE DRIVE CONTROLS

If driving with a traditional joystick is not the best option, alternative drive controls allow you to fully experience your chair's driving capabilities and benefits.



LINK-IT MOUNTING SYSTEM

Need to use additional switches or special controls? No problem, easily connect them using the universal Link-It Mounting System.

By using a combination of available modules and rods, you can create an endless number of combinations to accommodate what you need. Not only is it easy to position, Link-It also ensures the minimalistic, clean design of your wheelchair isn't compromised by additional controls. So whether it's an extra controller on the armrest, a push button on the backrest or a switch on the headrest, just link it.





QUICKIE



DESIGN AND APPEARANCE

DISTINGUISHED STYLE

With the Q500 series, we've taken style to the next level. With a matt-black chassis finish, hidden wiring and components combined with sharp-tailored upholstery, it's a sleek and contemporary design.

Fancy a change?

The Q500 series with SEDEO PRO, also come with a choice of 7 contemporary colour inlays to perfectly complement its striking matt-black bodywork. Simply replace the colour inlays. Styling your powered wheelchair has never been so simple.



TECHNICAL DATA / Q500 R - Q500 M - Q500 F WITH SEDEO PRO

	Q500 R	Q500 M	Q500 F		Q500 R	Q500 M	Q500 F
Max. user weight	160 kg	160 kg	160 kg	Max. speed	13 kph	13 kph	13 kph
Seat width	380-560 mm	380-560 mm	380-560 mm	Battery size	60 Ah; 80 Ah	60 Ah; 80 Ah	60 Ah; 80 Ah
Seat depth	400-600 mm	400-600 mm	400-600 mm	Electronics	VR2, R-Net	VR2, R-Net	VR2, R-Net
Seat height	415-490 mm	415-490 mm	415-490 mm	Turning radius	950 mm	560-630 mm	735 mm
Backrest height	440-620 mm	440-620 mm	440-620 mm	Max. reversing width	1258 mm	1100 mm	1240 mm
Total width	640 mm	610-620 mm	640 mm	Max. safe slope	10°	8°	8°
Total length	Max. 1180 mm	Max. 1280 mm	Max. 1340 mm	Max. kerb climbing	80 mm (100 mm with KC)	100 mm	70 mm
Total height	Max. 1400 mm	Max. 1400 mm	Max. 1400 mm	Drive wheels	14"	13"	14"
Powered recline	90°-150°	90°-150°	90°-150°	Castor wheel options	9"-10"	7"	9"-10"
Anti-shear back movmt.	100 mm (optional)	100 mm (optional)	100 mm (optional)	Colour inlay choices	7 Colour Accents	7 Colour Accents	7 Colour Accents
Tilt angle	30°/50°	30°/50°	30°/50°	Max. range (ISO 7176-4)	Max. 40 km	Max. 40 km	Max. 40 km
Seat functions	Tilt, Lift & Tilt, Recline	Tilt, Lift & Tilt, Recline	Tilt, Lift & Tilt, Recline	Wheelchair weight	From 133 kg	From 130 kg	From 133 kg
Crash tested	Yes - approved for forward facing transport in a vehicle (ISO 7176-19). Attachment: 4 chair fixing points with 2 double straps rear & 2 single straps front. Overall rating of belt restraint accommodation according ISO 7176-19:2008 Annex D=A. No anchored belt restraints allowed. The ease of access and manoeuvrability in vehicles depends on the wheelchair size.			Transport without user	Fix the wheelchair safe in the vehicle (e.g. car or plane).		
				Intended use and environment	For people with limited mobility; self-propelling or being pushed by an assistant, indoor and outdoor use. Wheelchair type class B.		



Sunrise Medical GmbH
 Kahlbachring 2-4.
 69254 Malsch/HD
 Tel.: +49 (0) 72 53 / 9 80-0
 Fax: +49 (0) 72 53 / 9 80-222
 E-Mail: international@sunrisemedical.de
www.SunriseMedical.eu

For further information on the full specification, options and accessories, please refer to the order form. All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.

If you are visually impaired, this document can be viewed in PDF format at SunriseMedical.eu.

