

GEMINO 30
carbon

GEMINO 30M
carbon



A rollator is a great way to maintain independence and live an active lifestyle. People who need a little support with balance and need to sit down often can benefit from a rollator. Rollators help you move around by offering physical support, safety, and stability. Perfect for seniors or younger people after injuries or with health conditions.



GEMINO 30 Carbon

The Premium lightweight Carbon Rollator for an Active Lifestyle

The Gemino 30 Carbon rollator is **very lightweight** with a **patented fold click system** and an elegant design. Thanks to the light weight, the rollator can be easily lifted for transportation or storage and makes it super maneuverable in tight spaces. This makes the rollator a perfect companion for everyday life and travel.

Carbon is extremely strong and solid yet lightweight, and therefore used wherever high demands are placed on the product. Despite the very low weight, the rollator is extremely stable and always provides a safe rolling experience indoors and outdoors – cobblestones or uneven paths are no problem.



// PATENTED FOLD-AND-CLICK SYSTEM

Thanks to the patented fold-and-click system the rollator folds to a compact size and automatically locks into place. Together with the low weight, the rollator can be easily lifted and stowed away – even with just one hand if necessary.



// WHEELS THAT MAKE THE DIFFERENCE

The large wheels enable walking comfortably – on any terrain.



// ANATOMICAL SHAPED HANDGRIPS

Ergonomically designed hand grips provide extra comfort and a secure feel for support while walking.

// HEIGHT ADJUSTABLE PUSH HANDLES

The simple height adjustment with memory function ensures an upright posture, minimizing a stooped posture.



// INCREASED SEATING SPACE

Wide parallel front tubes provide unobstructed space for comfortable seating.

// EASY HANDLING OF PURCHASES

A lot to carry? The advanced net bag can be easily opened with one hand, offers plenty of space for purchases and carries up to 5 kg.

The Gemini 30 Carbon is a loyal indoor and outdoor companion.

Tight spaces? With a width of only 600 mm, you can move around in narrow and small spaces. You can also safely maneuver in even tighter places by partially folding the rollator.

Need a break? The built-in seat combined with the back and wide space allows for a comfortable rest when needed. Lean back and relax.



PRACTICAL ACCESSORIES FOR YOU

Gemino rollators are not just easy to adjust, they come with a full range of attractive accessories to make everyday life easier. Choose the options that fit your needs: a stylish bag, a practical crutch holder, or a comfortable back support.

// COLOURFUL ADVANCED BAGS

Compliment your Gemino with a stylish, closable bag in your favourite colour. With handles and shoulder strap for easy carrying.



- attach & remove with one click. Simple one hand operation.



Black



Alcantara black



Graphite grey



Mint green



Red



Blue



Green



Pink



// BACK SUPPORT

Provides you support to sit comfortably and safely when it's time to take a break. No need to remove when folding the rollator.




// ADJUSTABLE COMFORT BACK SUPPORT

For additional comfort, choose this large, softly padded back support. Quick to adjust, easily folds with the rollator.



// TABLE TRAY

Practical tray with cup holder allows you to transport food and beverages safely and with ease.

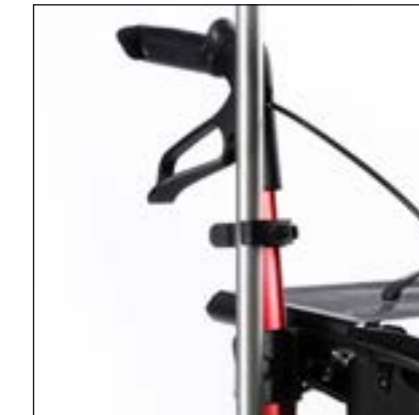


// CRUTCH HOLDER

Keeps your crutch or walking cane securely fastened when walking. Easy to mount, easy to access.



// UNIVERSAL BRACKET

Suitable for fitting your choice of personal accessories.

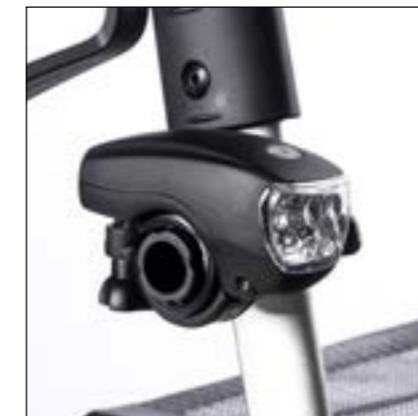


// LIGHT

Ensures you're well visible in the dark - for additional safety.



// ONE HAND BRAKE

Comfortable and safe usage with one hand.

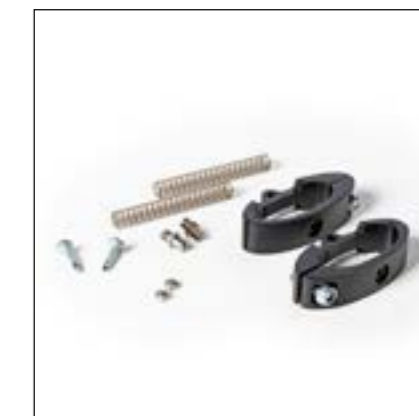

// REVERSE BRAKE KIT

The rollator immediately stops moving when the brake is released.

Technical data

GEMINO 30 Carbon // ROLLATOR

		Gemino 30 Carbon	Gemino 30M Carbon
Total height, from surface to top of pushing	mm	780 - 1000	730 - 880
Total length	mm	650	650
Total width	mm	600	600
Width between pushing handles	mm	470	470
Seat height	mm	620	550
Height folded	mm	800	760
Length folded	mm	650	650
Width folded	mm	230	230
Turning circle	mm	840	840
Front/rear wheels	mm	200 x 36	200 x 36
Total weight without basket	kg	5,6	5,4
Weight basket, standard	kg	0,3	0,3
Max. user weight	kg	150	130
Max. weight in basket	kg	5	5
Recommended user height	mm	1500 - 2000	1350 - 1700



Scan the QR code with your phone's camera to see more information.



Sunrise Medical GmbH
 Kahlbachring 2-4
 69254 Malsch/HD
 Tel.: +49 (0) 72 53 / 9 80-0
 Fax: +49 (0) 72 53 / 9 80-222
 E-Mail: export@sunrisemedical.de
www.SunriseMedical.eu

